



BIOPLAN[®]

FLEXIBLE
SYMPTOMATIC
FLATFOOT



CE 0426

 **BRM[®]**
Extremities

PRODUCT DESCRIPTION

Endorthesic implant for the flatfoot correction, through the mechanical action of filling the tarsal cavity, counteracting the sliding of the talus over the calcaneus and restoring the proper shape of the plantar arch.

Made of titanium alloy or PEEK with a titanium alloy marker for radiographic detection.

Available in 7 sizes.

ABSOLUTE CONTRAINDICATIONS

- Rigid flatfoot or with associated deformities
- Flatfoot with severe forefoot adduction
- Chronic ruptures of the posterior tibial tendon
- Symptomatic arthritis
- Evident joint degeneration, such as advanced subtalar joint osteoarthritis
- Overlying structural alignment deformities
- Peroneal muscle spasm
- Suspected or active infections

INTENDED USE

- Symptomatic flexible flatfoot in children and adolescents
 - In combination with calcaneal medialization osteotomies, in patients with significant hindfoot valgus deformity in the unstable pronated adult foot
- Excessive ligamentous laxity
 - Patient allergy to implant materials
 - Insufficient bone structures to ensure proper implant fixation
 - Local bone tumors
 - Weight exceeding manufacturer’s specified limit
 - Severe systemic diseases or metabolic disorders affecting the relevant anatomical area



For additional information and document downloads

IMPLANTS

	TITANIUM	Ø	L
	BRSBIOP0713V	7.0	13
	BRSBIOP0814V	8.0	14
	BRSBIOP0915V	9.0	15
	BRSBIOP1016V	10.0	16
	BRSBIOP1117V	11.0	17
	BRSBIOP1218V	12.0	18
	BRSBIOP1319V	13,0	19

	PEEK	Ø	L
	BRSBIOP0713P	7,0	13
	BRSBIOP0814P	8,0	14
	BRSBIOP0915P	9,0	15
	BRSBIOP1016P	10,0	16
	BRSBIOP1117P	11,0	17
	BRSBIOP1218P	12,0	18
	BRSBIOP1319P	13,0	19

INSTRUMENTATION



UAOI0MA000 Quick Coupling Handle



UAOI0SA00S AO quick-coupling silicone handle



UAOI0SAC0S AO quick-coupling silicone ratchet handle



SBPI1CI001 Inserter Screwdriver



SBPI1CE002 Extractor Screwdriver



SBPI10S001 Probe



SBPI1SF001 Push wire



SBPI10P007 Trial size 7



SBPI10P008 Trial size 8



SBPI10P009 Trial size 9



SBPI10P010 Trial size 10



SBPI10P011 Trial size 11



SBPI10P012 Trial size 12



SBPI10P013 Trial size 13




SBPI1PG001 Guide wire Ø 2 mm



SBPI000000B Box BioPlan

KIT BIOPLAN

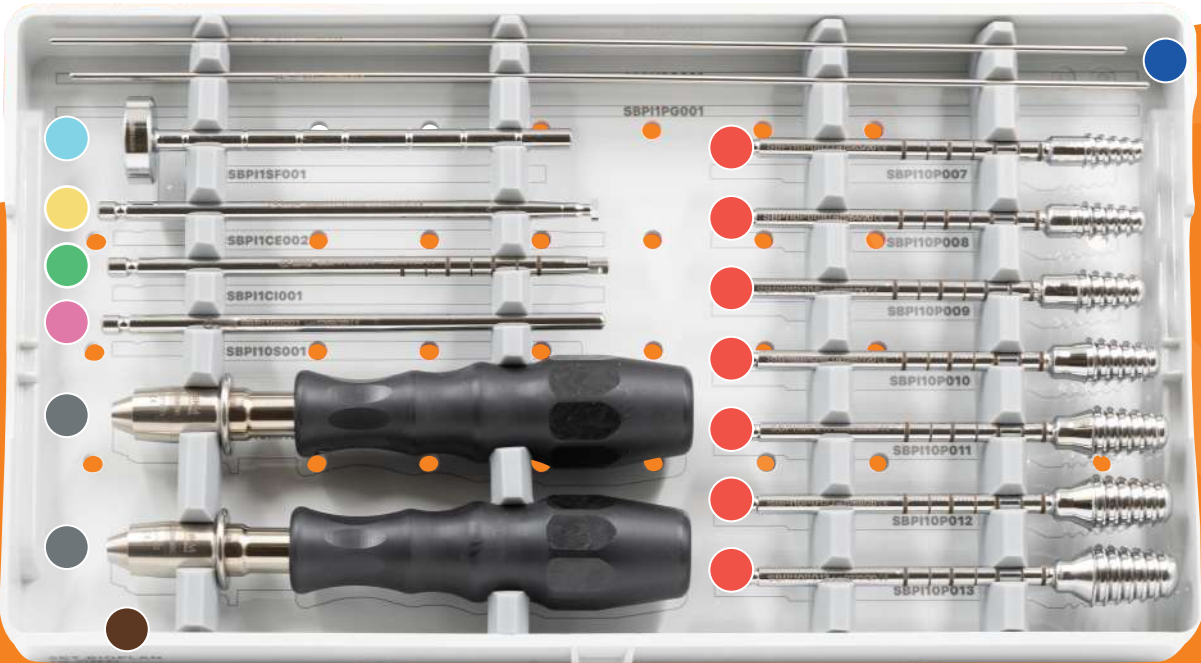

Distributed by:



BRM EXTREMITIES Srl

Via Lorenzo Mascheroni 29
20145 Milano (MI) Italy
T +39 02 89415508
brm-extremities.com

Headquarter

Via Privata ai Campi 7
23862 Civate (LC) Italy
T +39 0341 1693087
info@brm-extremities.com

D-03A-001-A-EN-REV00