

BRM TOOL VIP SCREWS

VERIFIED
INTRATHREAD
PERFORMANCE




CE 0426

PRODUCT DESCRIPTION

Interfragmentary stabilization screws with chamfered proximal thread for osteosynthesis of small bone fragments. Combine great self-tapping and self-drilling properties with great fracture reduction capabilities, to provide optimal internal fixation.

Available in 3 different diameters and different thread lengths.

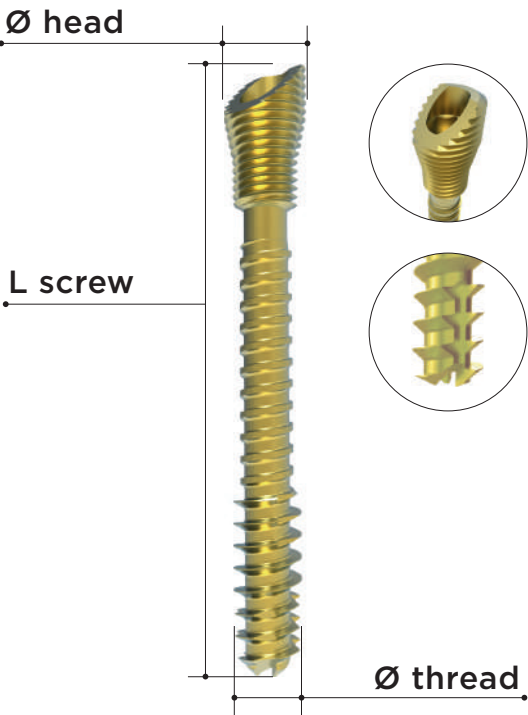

CLINICAL INDICATION

Fixation of small and medium sized bone fragments:

- epiphyseal detachments
- osteotomy for the correction of hallux valgus
- distal osteotomy of the second, third or fourth metatarsal for painful plantar calluses
- ankylosis of the ankle
- ankylosis of the wrist
- metacarpal and phalangeal fractures

45° CHAMFERED HEAD
SLIGHT INTRATHREAD
SELF-DRILLING
AND SELF-TAPPING
TITANIUM ALLOY Ti6AL4V



For additional information and document downloads

Diameter	2.5	3.0	4.0
L screw	From 16 to 40 mm	From 16 to 60 mm	
Connection	Hex 1.7	Hex 2.0	Hex 2.5
K-Wire	Ø 1.0		Ø 1.5
Thread	Single progressive		
Color code	Green	Yellow	Blue

STERILE SCREWS 2.5



Reference	Ø	L
SSL125016S	2.5	16
SSL125018S	2.5	18
SSL125020S	2.5	20
SSL125022S	2.5	22
SSL125024S	2.5	24
SSL125026S	2.5	26
SSL125028S	2.5	28
SSL125030S	2.5	30
SSL125032S	2.5	32
SSL125034S	2.5	34
SSL125036S	2.5	36
SSL125038S	2.5	38
SSL125040S	2.5	40

STERILE SCREWS 3.0



Reference	Ø	L
SSL130016S	3.0	16
SSL130018S	3.0	18
SSL130020S	3.0	20
SSL130022S	3.0	22
SSL130024S	3.0	24
SSL130026S	3.0	26
SSL130028S	3.0	28
SSL130030S	3.0	30
SSL130032S	3.0	32
SSL130034S	3.0	34
SSL130036S	3.0	36
SSL130038S	3.0	38
SSL130040S	3.0	40
SSL130042S	3.0	42
SSL130044S	3.0	44
SSL130046S	3.0	46
SSL130048S	3.0	48
SSL130050S	3.0	50
SSL130052S	3.0	52
SSL130054S	3.0	54
SSL130056S	3.0	56
SSL130058S	3.0	58
SSL130060S	3.0	60

TOOL VIP
SCREWS

STERILE SCREWS 4.0



Reference	Ø	L
SSL140016S	4.0	16
SSL140018S	4.0	18
SSL140020S	4.0	20
SSL140022S	4.0	22
SSL140024S	4.0	24
SSL140026S	4.0	26
SSL140028S	4.0	28
SSL140030S	4.0	30
SSL140032S	4.0	32
SSL140034S	4.0	34
SSL140036S	4.0	36
SSL140038S	4.0	38
SSL140040S	4.0	40
SSL140042S	4.0	42
SSL140044S	4.0	44
SSL140046S	4.0	46
SSL140048S	4.0	48
SSL140050S	4.0	50
SSL140052S	4.0	52
SSL140054S	4.0	54
SSL140056S	4.0	56
SSL140058S	4.0	58
SSL140060S	4.0	60




TOOL VIP
SCREWS

INSTRUMENTATION



- UAOI0MA000** Quick Coupling Handle
- UAOI0SA00S** AO Quick Coupling Silicone Handle

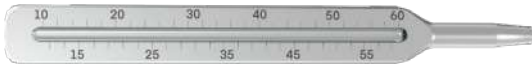

- STOI232000M** Depth gauge 60 mm



- STVI220XL0P** Long cannulated drill bit Ø 1.8 mm for screws Ø 2.5 mm
- STVI224XL0P** Long cannulated drill bit Ø 2.0 mm for screws Ø 3.0 mm
- STVI234XL0P** Long cannulated drill bit Ø 2.9 mm for screws Ø 4.0 mm



- STVI2250BK** Basic countersink Ø 2.5 mm for screws Ø 2.5 mm
- STVI2300BK** Basic countersink Ø 3.0 mm for screws Ø 3.0 mm
- STVI2400BK** Basic countersink Ø 4.0 mm for screws Ø 4.0 mm



- STVI2ES017S** Slim cannulated screwdriver Hex 1.7 for screws Ø 2.5 mm
- STVI2ES020S** Slim cannulated screwdriver Hex 2.0 for screws Ø 3.0 mm
- STVI2ES025S** Slim cannulated screwdriver Hex 2.5 for screws Ø 4.0 mm



- STVI2ES017P** Solid screwdriver Hex 1.7 for screws Ø 2.5 mm
- STVI2ES020P** Solid screwdriver Hex 2.0 for screws Ø 3.0 mm
- STVI2ES025P** Solid screwdriver Hex 2.5 for screws Ø 4.0 mm



- UKW100150T** Tube for K-wire sterilization


- STVIEG000BP** Sterilization tray for instruments Tool Vip

K-WIRES



- STOIKW10150803N** K-Wire Ø 1.0 x 150 mm with marker at 80 mm Pack of 3 pcs - Non Sterile for screws Ø 2.5 - 3.0 mm
- STOIKW10150810N** K-Wire Ø 1.0 x 150 mm with marker at 80 mm Pack of 10 pcs - Non Sterile for screws Ø 2.5 - 3.0 mm
- STOIKW15150803N** K-Wire Ø 1.5 x 150 mm with marker at 80 mm Pack of 3 pcs - Non Sterile for screws Ø 4.0 mm
- STOIKW15150810N** K-Wire Ø 1.5 x 150 mm with marker at 80 mm Pack of 10 pcs - Non Sterile for screws Ø 4.0 mm

KIT TOOL VIP

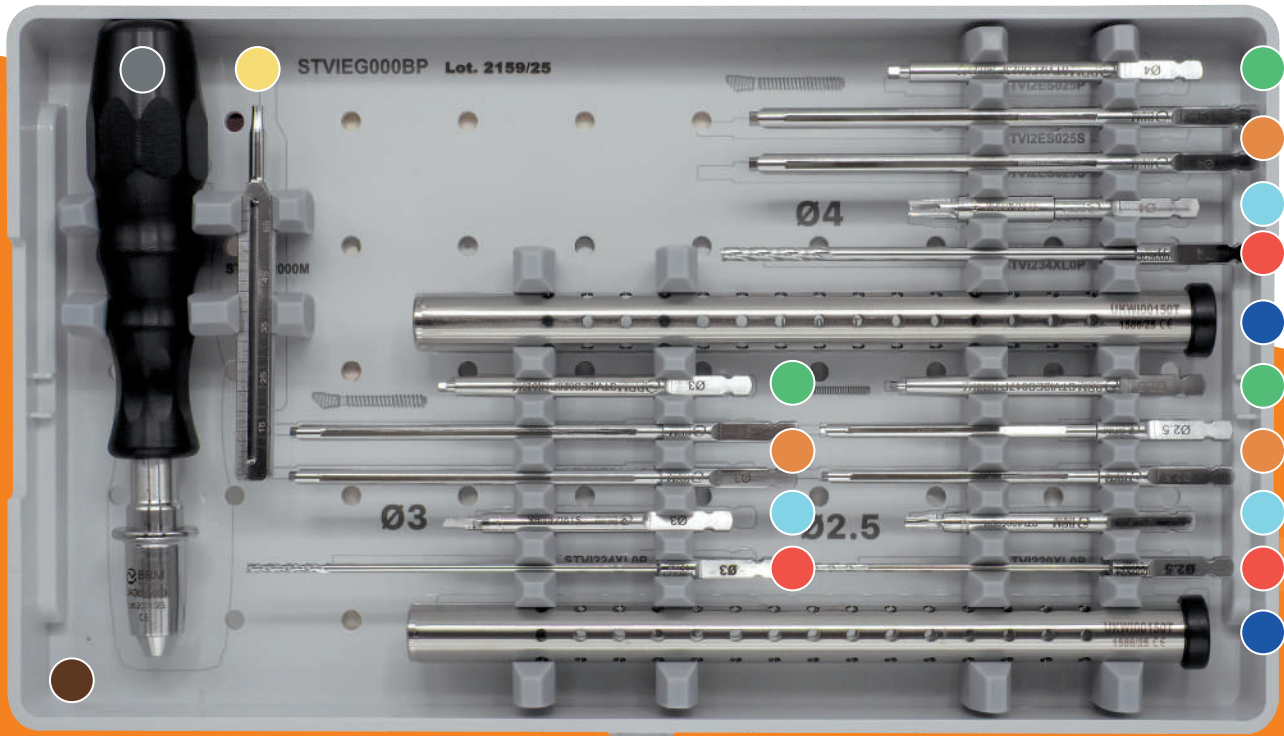

Distributed by:



BRM EXTREMITIES Srl

Via Lorenzo Mascheroni 29
20145 Milano (MI) Italy
T +39 02 89415508
brm-extremities.com

Headquarter

Via Privata ai Campi 7
23862 Civate (LC) Italy
T +39 0341 1693087
info@brm-extremities.com

D-01A-004-A-EN-REV00